


Order form - select the correct Cover


Flex **Neo Max BK**

XT
COVER

Prosthetic socket circumference at the very end of the cover in [mm]

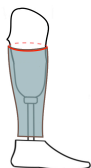


To take the measurement correctly, please mark the height of the cover in several places and place a ruler below the marking. The measurement should be made very carefully because moving the ruler even by a few millimetres up or down causes an incorrect reading of the circumference even by several centimetres.

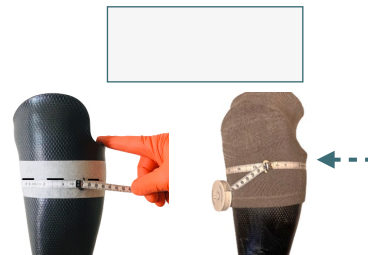


Height in [mm]

We measure from the level of the foot shell to the top of the knee. The gap between the cover and the foot will be taken by us, depending on the type of foot.

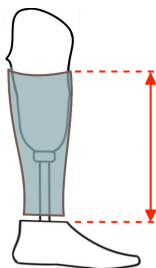

1

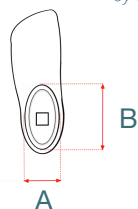

A [mm]


B [mm]

Foot manufacturer, type and size

For example : OB Taleo 1C50 L29

The silicone sleeve - option

Cover below the silicone sleeve or without

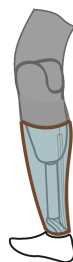

Cover over the silicone sleeve i.e. the sleeve is between the socket and the cover



Adapter

30

34

e.g.:
4WR95=3

Matte and semi-matte colours - tattoos unavailable

- | | | | |
|-----------------------------|-----------------------------|-----------------------------|-----------------------------|
| <input type="checkbox"/> 21 | <input type="checkbox"/> 41 | <input type="checkbox"/> 61 | <input type="checkbox"/> 81 |
| <input type="checkbox"/> 23 | <input type="checkbox"/> 43 | <input type="checkbox"/> 63 | <input type="checkbox"/> 83 |
| <input type="checkbox"/> 25 | <input type="checkbox"/> 45 | <input type="checkbox"/> 65 | <input type="checkbox"/> 85 |

Shiny colours & matt varnish

- | | | |
|-------------------------------------|---|-------------------------------------|
| <input type="checkbox"/> Skin Tone1 | <input type="checkbox"/> Skin Tone2 | <input type="checkbox"/> Skin Tone3 |
| <input type="checkbox"/> Graphite | <input type="checkbox"/> Grey | <input type="checkbox"/> White |
| <input type="checkbox"/> Black | <input type="checkbox"/> Matt varnish (charged extra) | |

Tattoo (charged extra)

Available only for Skin Tone3, White and Grey colour

Catalog number e.g.: 08-1001

Own

Check list :


Photo of the foot or shell.

Both the outer and inner edges of the shell contour must be visible. The camera lens must be perpendicular to the pyramid.



4 Photos The camera must be placed perpendicular to the legs, in a low position, and both legs must be centered on the photo. Please remember that all prosthetic components must be completely visible, including the socket. Please remove shoes, socks and any object that covers the prosthesis.


Patients do not need to remove the socks from the healthy leg.

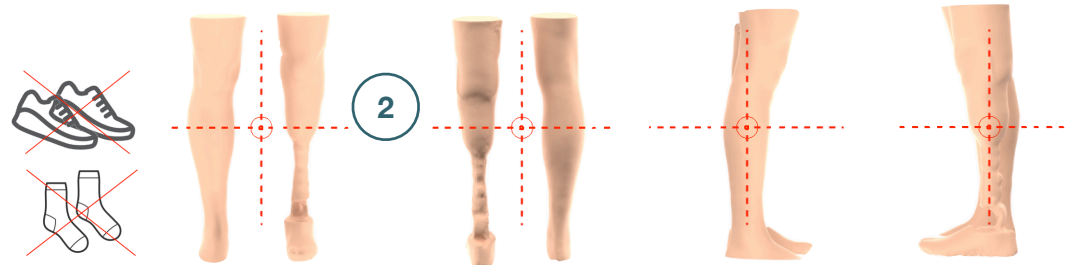

Scan of the prosthesis and the patient's leg



Formats
.stl
.obj


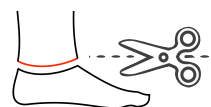


- All our covers can be shortened from the bottom by several mm.
- In the factory state, the lower part of the cover is flat.
- In the case of the foot as shown in the picture, the lower part of the cover can be ground to the shape of the foot.


Send the form and files to nextt@nextt.pl or use an external service such as **WeTransfer**

VERY IMPORTANT information - do not send.

1

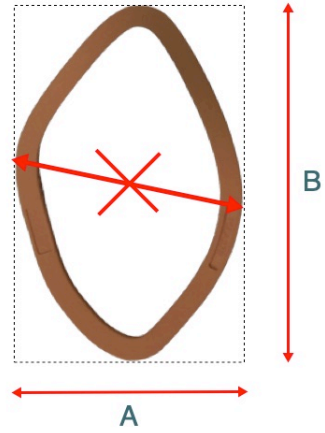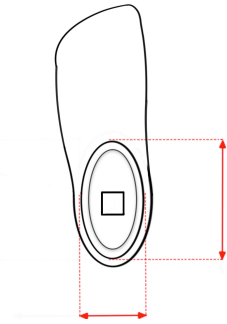

2

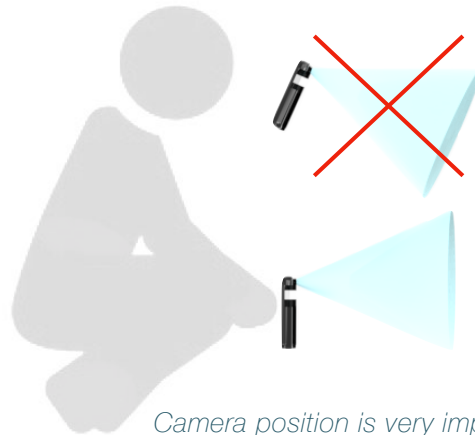

Camera position is very important

COLOURS - do not send.

Matte and semi-matte colours - tattoos unavailable



Shiny colours - tattoos available

