


# PHANTOM

ABOVE KNEE AMPUTATION, TRANSFEMORAL

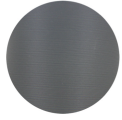

***Phantom***

## TECHNICAL DETAILS

- Available for all knee joints.
- Waterproof.
- Size : Any.
- Covers can be shortened by a dozen mm.
- Compatible with the knee rotation adapter.
- Shape : anatomical, made according to a scan of a healthy leg.
- Material : soft and flexible, highly resistant to water and bacteria.
- Adapted to kneeling, if the knee prosthesis design allows for it.
- The top cover is attached to the prosthetic socket with a special button, allowing you to remove the cover at any time.
- The bottom cover has magnets, allowing you to put on and take off the cover quickly.
- The cross-section of the bottom part of the bottom cover is similar to the shape of the foot shell.
- Painless Tattoos - extra paid option.
- Matte version - extra paid option.
- There is a gap of approximately 7 mm between the bottom cover and the foot. For very soft feet, such as carbon fiber feet, the distance should be significantly greater than 7 mm.




## COMPARISON

**Duo**

**Phantom**

KNEE SHAPE

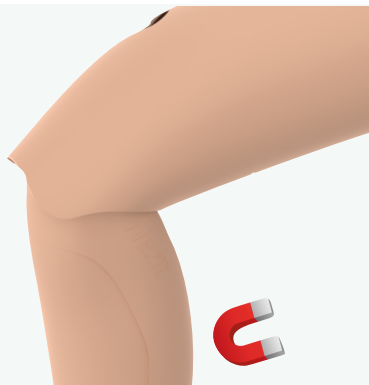

BOTTOM COVER

**Flex Neo AK**

SIMILAR TO

**MAX AK**

TOP COVER ATTACHMENT METHOD

PERMANENTLY USING DOUBLE-SIDED TAPE

EASY TO REMOVE THANKS TO OUR PROPRIETARY BUTTON

ADAPTATION TO KNEELING

YES, AS LONG AS THE DESIGN OF THE KNEE ALLOWS FOR IT

MATERIAL

SOFT AND FLEXIBLE, HIGHLY RESISTANT TO WATER AND BACTERIA

# PRODUCT DESCRIPTION

The Phantom Cover is our latest achievement in flexible cosmetic covers. It has been designed to closely mimic a healthy limb while walking and, most importantly, sitting down when the prosthesis is bent.

Our proprietary solution uses a unique "button" and magnets, making it easy to put the covers on and take them off (both the bottom and the top ones).

easy to put on clothes over the covers.

Compatible with the knee rotation adapter.

Highly resistant to mechanical impact, they effectively protect the prosthesis. It is especially recommended for very active people. In addition, the material from which they are made is very durable, resistant to water and both low and high temperatures.




Adapted to kneeling





Adapted to sitting


# COLOURS

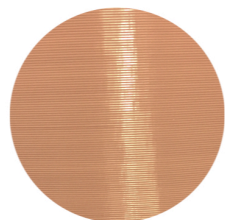


Skin tone

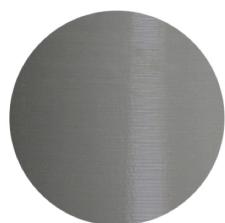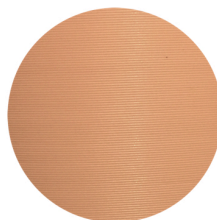



Grey

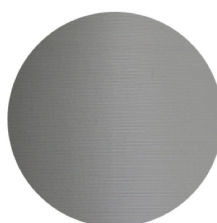


White

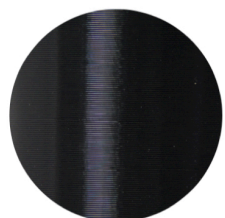

Black

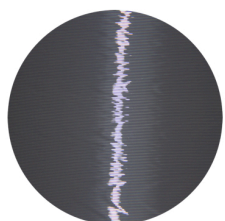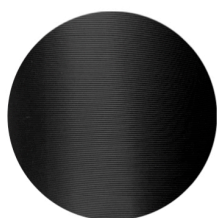

Graphite

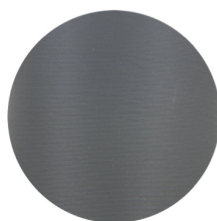

Glossy

Matte